
Sonno

Modello matrimoniale
Model double

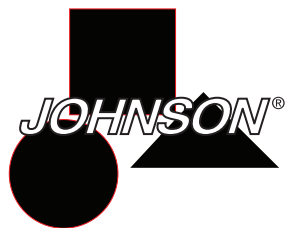Modello singolo
Model single

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GB GENERAL WARNINGS

Please read the instructions for use and general warnings carefully before using this appliance. Keep this manual for future reference. The following safety and danger indications not only serve to protect the health of those using the appliance and that of third parties, but also to protect the product.

CAUTION:

This bed warmer must not be used by persons insensitive to heat or other highly vulnerable people unable to react to overheating.

Not suitable for use with children under 3 years of age. Children under 3 years of age must not use this bed warmer due to their inability to react to overheating. This product must be kept out of reach of children under 3 years of age.

This bed warmer must not be used for small children (3-8 years of age), unless the switch has been set by an adult and the child has been sufficiently instructed on the safe use of the bed warmer.

The bed warmer may be used by children of more than 8 years of age, or by people with limited physical, sensory or mental abilities, or by those without experience or necessary knowledge, provided they are under the super-

vision of adults or after they have received instructions regarding the safe use of the device and the understanding of the dangers inherent to it. Children must not play with the bed warmer. Cleaning and maintenance of the bed warmer to be carried out by the user must not be carried out by children without supervision. The appliance and its power cord must be kept out of reach of children under 8 years of age. Persons with pacemakers should not operate this device when they are in bed: it can be used to preheat the bed, but must be switched off and disconnected from the power supply before going to bed.

Consult your doctor before use.

Before connecting the appliance, check that the plate data correspond to those of the electrical distribution network. In case of incompatibility between the socket and the plug, contact qualified personnel for the necessary adjustment operations.

Do not operate the temperature control with with wet or damp hands and feet.

Always disconnect the power supply before cleaning or maintenance and if the appliance is not in use.

If the appliance's power cord is damaged, return the entire appliance to an authorized tech-

nical assistance center to replace the control.

If the cord cannot be removed, the bed warmer must be disposed of.

Recommendation: do not leave the bed warmer in use unnecessarily.

Do not place the the control under the pillow or between sheets or blankets when the appliance is connected.

Check the appliance frequently for signs of wear or deterioration. If there are signs of wear or deterioration, or if the appliance is damaged or has obvious malfunctions, contact an authorized technical assistance centre before use.

The appliance can be used to heat the bed, but must be switched off and disconnected from the power supply before going to bed.

CAUTION

Place the bed warmer only between the mattress and the first sheet and not over the persons in bed.

If the appliance is used on an adjustable bed, check that the bed warmer and the power cord are not trapped or blocked, for example in the hinges.

Do not operate the bed warmer when there are obvious folds.

Do not insert pins or needles. Do not use pins,

clothespins or other metal tools to hold the blanket in place, do not scratch the surfaces in any way.

The bed warmer cannot be used for medical purposes.

If the thermal protection device of the product activates, the appliance must be taken to an authorized technical assistance centre.

Do not wash or dry clean. Do not iron. Clean using a damp cloth only. Make sure that the blanket is completely dry before use.

To avoid damaging the insulating layer of the heating cables inside the product, it is absolutely forbidden to exposure to sunlight, rain or extreme temperatures.

Do not use in combination with other heat sources.

The product must be put into operation only via its regulator, without the aid of other devices, such as timers for automatic power-on. The product is not indicated for use in hospitals.

Do not use with water bed mattresses.

Store the bed warmer in a dry place.

When storing the appliance, allow it to cool before folding.

Do not place objects on the appliance during storage, which would cause creases and/or de-

form the appliance.

Do not place heavy objects (e.g. suitcases or laundry baskets) on the appliance, but only those with coverage function.

Do not use the bed warmer to keep pets warm.

Do not use foled.

In case of failure or malfunction of the appliance, switch it off and do not tamper with it. For any repairs, contact only an authorized technical assistance center.

Do not pull the power cord or the appliance itself to remove the plug from the socket.

ATTENTION:

Warranty does not cover any plastic parts of the product.

Warranty does not cover any glass parts of the product.

ATTENTION:

The damages of the power supply cord, originate by the wear aren't covered by guarantee; the repair will be in charge to the owner.

ATTENTION

If the need should arise to take or send the appliance to an authorised service centre, make sure to clean all its parts thoroughly. For hygienic purposes, the service centre shall reject appliances that are not perfectly clean on the outside or inside without making any controls or repairs.

The symbol  on the product or on its packaging indicates that this product may not be treated as household waste. Instead it shall be handed over to the applicable collection point for the recycling of electrical and electronic equipment.

By ensuring this product is disposed of correctly, you will help prevent potential negative consequences for the environment and human health, which could otherwise be caused by inappropriate waste handling of this product.

For more detailed information about recycling of this product, please contact your local city office, your household waste disposal service or the shop where you

purchased the product.

COMPONENTS

- A. Command group
- B. Temperature selector
- 0. Electric blanket OFF
- 1. First heat level
- 2. Second heat level

INSTRUCTIONS FOR USE

Positioning of electric blanket.

- The electric blanket should be positioned on a fixed and resistant mattress

- Position the underblanket with the temperature selection connector in contact with the mattress.

Take care the blanket is well spread over the mattress in the zone not covered by pillows. Completely unwind the connector cable and power supply cable.

- The temperature selector (A) of the heating element must be close to the shoulders for ease of use, so as to correctly use the heating adjustment function.

. (Folded parts cause high temperature so they can become dangerous).

- The socket of the heating element (D) should be near your shoulders for easy use and to correctly use the heating function.

. (Folded parts can cause a high temperature which can be dangerous).

- To always ensure the electric blanket is spread, clean and fold-free, always position the sheet too.

(Do not place other blankets or covers in between, as they could obstruct heat diffusion).

- The temperature selector (A) should be left outside, and not under the pillow, or any other object to avoid exerting pressure on it.

The electric blanket should be heated 30 minutes to 1 hour in advance before to lie down.

Ensure the switch (B) is on 0. Insert the connection plug in a suitable electric socket. Switch on the appliance by bringing the switch (B) to 1 or 2 according to the temperature desired.

The appliance is set to work 1 hour for maximum power.

TEMPERATURE SELECTION

Exert pressure upwards to select the two different temperatures:

- 0 Off-
- 1 First heat level
- 2 Second heat level.

SAFETY PRECAUTIONS

Before use, check the inner heating cable is not faulty or twisted, if so, do not use the appliance and repair it in a support centre.

Often check the status of the appliance. If you detect one of the conditions listed below, stop use and dispose of it.

- Power cable insulation stripped or damaged.

- Internal resistor sticking out of electric blanket.

SPECIAL WARNINGS

Place the underblanket only between the mattress and the sheet and never over the person in the bed.

The underblanket is not suitable for use in hospitals and not recommended on beds with adjustable mechanisms and joints that could damage the cable. Avoid creating creases and store it, even when not in use, in a dry and clean environment.

Do not insert pins, needles or other foreign bodies into the appliance's tissues. Do not dry clean or immerse in liquids. Do not iron or put in contact with hot elements. Do not lift or move the underblanket using the power cord.

- Do not use for children, people who are unable, insensitive to heat, or who have drunk alcohol that causes sleep.

- Do not use in conjunction with other heating appliances

- Use only as a lower deck, before each use it must be completely unfolded.

- Do not insert pins or anything else to position the blanket, this can damage the appliance or cause injury.

- Do not use if damp or wet.

Note: remove the plug from the socket without pulling the power cord or the connection cable.

- Do not wash it with water or dry

- Do not tamper with, repair or modify the appliance, do not replace the power cable. In case of problems, contact the manufacturer or an authorized service center.

STORAGE

When the electric blanket is no longer used, keep it as follows:

when storing the appliance, allow it to cool before folding it

wait for the electric blanket to be completely dry if previously cleaned

fold the electric blanket like a normal blanket

for a better storage it is advisable to store the electric blanket in the package supplied with the purchase

store the electric blanket in a dry place

Attention: do not place heavy objects on top of the electric blanket, when it is stored, so as not to damage it

CLEANING

To clean the temperature selector use a dry and soft cloth. This should be done always with the plug disconnected from the network. Do not wet or immerse the temperature selector and the cables in water or other liquids. To clean the heating element follow scrupulously the instructions on the product labels.

Any cleaning (always with the plug unplugged) should be made with very gentle suction and then passed with a sponge or damp cloth.

Before replacing the appliance for the use wait until it is completely dry.

Do not wet or immerse the appliance in water absolutely.

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GUARANTEE

The duration of the warranty runs from the date of purchase of the appliance, proven by the retailers stamp, and terminates after the period specified, even if the appliance was not used. The warranty covers all replacements or repairs necessary due to material or manufacturing defect. The warranty expires when the appliance is tampered with or when the defect is due to free of charge for goods delivered to our support centres. The client must pay the transport expenses only. The parts or appliances replaced become our property.

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